
Simple concentration monitoring system. Never miss a concentration change with real-time readings!



# Digital Immersion Refractometer PAN-1



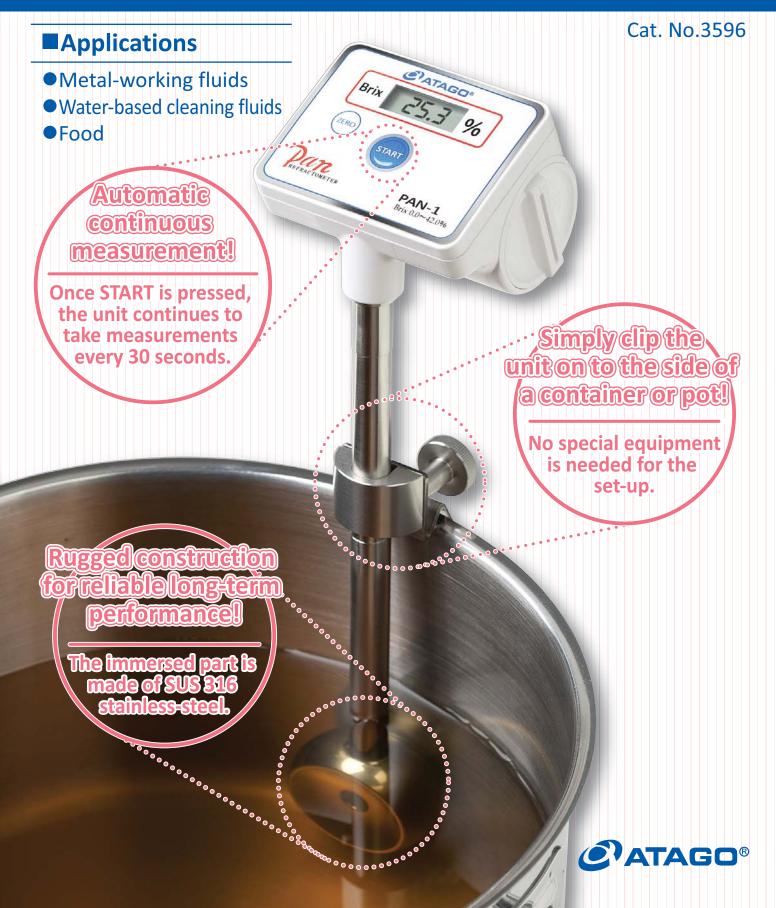

values are large and

easy to read.

ZERO Key

Immerse the prism

this to zero-set.

Press this once,

and continuous

measurement starts.

% Use water at a similar temperature to the liquid to be measured.

80mm

96

ABS

SUS316L

Optical glass

SUS316L

Silicon(packing)

START Key

(Immersed part

Prism

Dimensions

300mm

Battery Case

LSD NiMH.

This can attach to

various shapes of


containers.

Clamp

Size D battery.

Compatible with

## Measurement Method

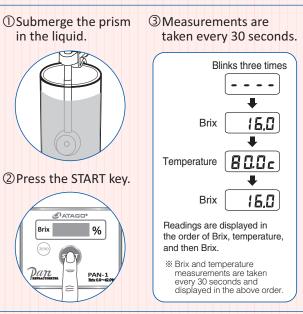

## Displayed Value

This instrument measures the refractive index of a liquid and displays the converted value in "Brix," which is a percentage measurement of sugar in water. The concentration of a substance other than sugar can be calculated by creating a conversion table of the desired concentration % and corresponding Brix %. Contact ATAGO for further details.

## Specifications

PAN-1	Cat. No.3596
Measurement range	Brix 0.0 to 42.0% Temperature 10.0 to 99.9°C
Resolution	Brix 0.1% Temperature 0.1°C
Measurement accuracy	Brix ±0.2% Temperature ±0.5°C
Measurement temperature	10 to 95°C
International Protection class	Immersed part IP67 Display IP65
Dimensions and weight	80(W)×300(L)×72(H)mm 610g (main unit only)
Battery life	2 months(when an alkaline battery is used)
Power supply	Size D alkaline battery × 1 Compatible with LSD NiMH.
Ambient temperature	10 to 45°C
Accessories	<ul><li>①Size D battery</li><li>②Instruction manual</li></ul>

#### All ATAGO refractometers are designed and manufactured in Japan.

### ATAGO CO., LTD.

Headquarters: The Front Tower Shiba Koen, 23rd Floor 2-6-3 Shiba-koen, Minato-ku, Tokyo 105-0011, Japan TEL : 81-3-3431-1943 FAX : 81-3-3431-1945 overseas@atago.net http://www.atago.net/

CATAGO U.S.A., Inc. TEL : 1-425-637-2107 customerservice@atago-usa.com

CATAGO INDIA Instruments Pvt.Ltd. TEL: 91-22-2833-8038 / 8076 customerservice@atago-india.com

TEL : 662-982-8718-9 customerservice@atago-thailand.com

CATAGO BRASIL Ltda. TEL : 55 16 3916-6000 customerservice@atago-brasil.com

ATAGO ITALIA s.r.l. TEL: 39 2 36557267 customerservice@atago-italia.com CATAGO CHINA Guangzhou Co., Ltd. TEL : 86 20-38106065 info@atago-china.com



Fastening part

Adjustable height range

69mm

152mm

### \* Specifications and appearance are subject to change without notice.

HACCPIGMPIGLP

ATAGO products comply with HACCP, GMP, and GLP system standards.